
medisana®

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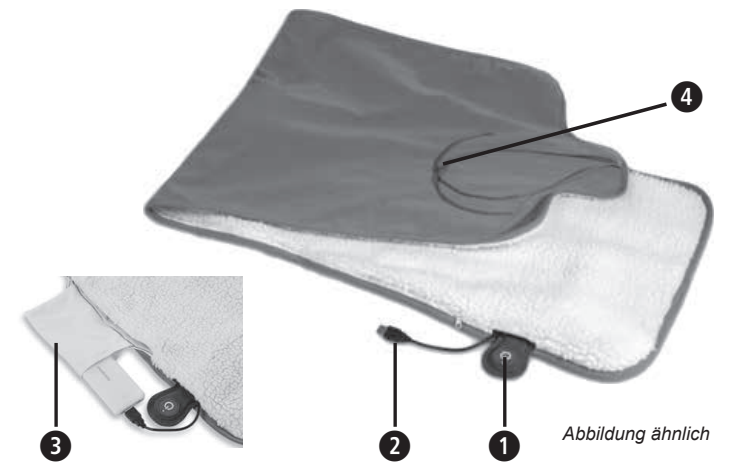
medisana®

GB Instruction manual
Heated Pad OL 800 / OL 850

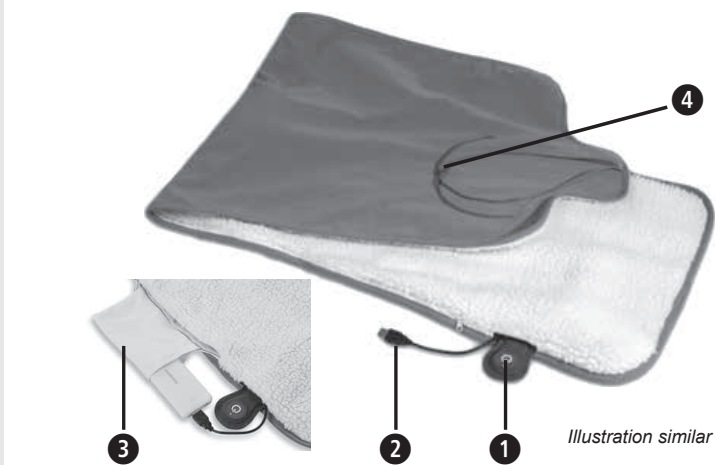- ❶ On/Off button
- ❷ USB plug for power bank*
- ❸ Integrated pocket for power bank*
- ❹ Strap to roll up

* Powerbank not included in scope of delivery

Explanation of symbols

- Do not insert needles into the device!
- Not to be used by very young children (0-3 years)!
- The cover can be washed by hand! Follow the instructions on the label!
- Do not bleach!
- Do not dry the device in a tumble drier!
- Do not iron!
- Do not wring!
- Do not dry clean

IMPORTANT!
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes provide useful additional information on the installation or operation.

LOT LOT number

Manufacturer

Keep dry

Recycling symbols/codes: These are used to provide information about the material and its proper use and recycling.

GB Safety instructions
 Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Power supply

- Before connecting the unit to the power supply, make sure that it is switched off and that the specified electrical voltage is observed.
- Always switch off the appliance after use with the button on the control unit and disconnect the power supply.
- Never carry, pull or turn the device by the cable.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.

For special people

- Do not use the device on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Children under the age of three are not allowed to use this appliance due to their inability to react to overheating.
- The appliance is not allowed to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or unless the child has been adequately instructed on how to operate the controls safely.
- If you have any health concerns, talk to your doctor before use.

Operating the device

- Only use the device for its intended use as in the instruction manual.
- If misused, all guarantee rights shall become null and void.
- Before each use, carefully check the article for signs of wear and/or damage. Do not operate the appliance if you notice wear, damage or signs of improper use.
- If an error occurs, check that the power supply is properly established. Do not open the unit. There are no user-serviceable parts inside the unit.
- The Heated Pad must not be operated folded or bent.
- Do not immerse the device in water or other liquids. If liquid does enter the device, disconnect the power supply immediately.
- Do not use the device in a wet environment (e.g. in the bathroom or similar).
- Do not keep or place any electrical device on a surface where it is in danger of falling into a bathtub or sink.
- Do not use the device while you are sleeping or lying in bed. Danger of burns or heat stroke!
- The device is not suitable for medical use (e.g. in hospitals).
- Very long periods of use can cause overheating. Always allow the appliance to cool down before using it again.
- Never operate the appliance unattended.
- Do not attach safety pins or other pointed or sharp objects to the appliance or prick it.

- Do not cover the unit when it is on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury.
- Never store or use the appliance near an electric cooker or other heat sources.

Maintenance and cleaning

- Do not repair the unit yourself in the event of a fault. Repairs may only be carried out by an authorised dealer or a suitably qualified person.
- Cleaning and user maintenance shall not be made by children without supervision.
- When storing the appliance, allow it to cool down before storing it.
- Do not place any objects on top of the item during storage to avoid sharp kinks.
- Do not crease the appliance by placing items on top of it during storage.
- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, contact the customer service before switching it on again.

Scope of supply

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the it and contact your dealer or your service centre.

The following parts are included:

- 1 **medisana** Heated Pad **OL 800 / OL 850**
- 1 Storage pouch
- 1 Short instruction manual

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Operating

The **medisana** Heated Pad **OL 800 / OL 850** is characterised by wireless use thanks to the power bank (not included). For operation, you need a commercially available power bank with a USB socket with an output power of 5V/2A (e.g. the **medisana** power bank art. 60279).

- Connect the USB plug **❷** on the device to the USB socket on the power bank.
- Lay out the Heated Pad according to your needs.
- Press the on/off button **❶** to switch the unit on. The heat switches on at the lowest level and the LED lights up in yellow.
- By pressing the on/off button **❶** for a long time (hold for approx. 1 second), you can switch between heating levels 1, 2 or 3 (medium level 2 = LED orange; highest level 3 = LED red).
- The appliance switches off automatically after 20 minutes or can be switched off by short pressing the the on/off button **❶** at any time.

Stow the powerbank in the integrated pocket **❸**.

Power bank (not included in scope of delivery)

The power bank to be used must be fully charged before first use.

Troubleshooting

If there is no function AT ALL after switching on, please check whether the power supply is correctly established. If this is the case, the device is defective. Please do not use the device and contact the service center.

If the LED flashes, this is an indication of other possible faults:

- LED flashes 3 times - short circuit in the heating wire.
- LED flashes 7 times - circuit not closed
- LED flashes 5 times - short-circuit or interruption of the circuit NTC sensor

If such errors occur, disconnect the warming blanket from the power supply, wait 5 minutes and try to start it up again. If the error cannot be eliminated, stop using the blanket and contact the service centre.

Cleaning and maintenance

- Before cleaning, disconnect the power supply and allow the appliance to cool down for at least 10 minutes.
- Hand wash: Detach the removable cover and soak it in mild soap and cold water for about 15 minutes. Then rinse in fresh, cold water.
- Keep the cable and control unit away from heat, hot surfaces, moisture and liquids.
- Do not wring or rub the product during or after cleaning!
- May not be bleached, ironed or dry cleaned!
- Do not use the article again until it is completely dry.
- Untwist the cable if it is twisted.
- For storage, the item can be rolled up and tied with the cord **❹** and packed in the storage bag. Make sure that the item has cooled down completely before storing it.

Technical specifications

Name and model: **medisana** Heated Pad **OL 800 / OL 850**

Automatic switch-off: 20 min.

Input: 5V== 10W

Storage conditions: Cooled down, clean and dry

Dimensions approx.: 135 x 45 cm

Weight approx.: 0.6 kg

Art. number: 60281 (OL 800)
60282 (OL 850)

EAN number: 40 15588 60281 8
40 15588 60282 5

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH, Carl-Schurz-Str. 2,
41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/60282>.

Information about service can be found here:
www.medisana.com/servicepartners

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